


Rocket® Seldinger Chest Drainage Set

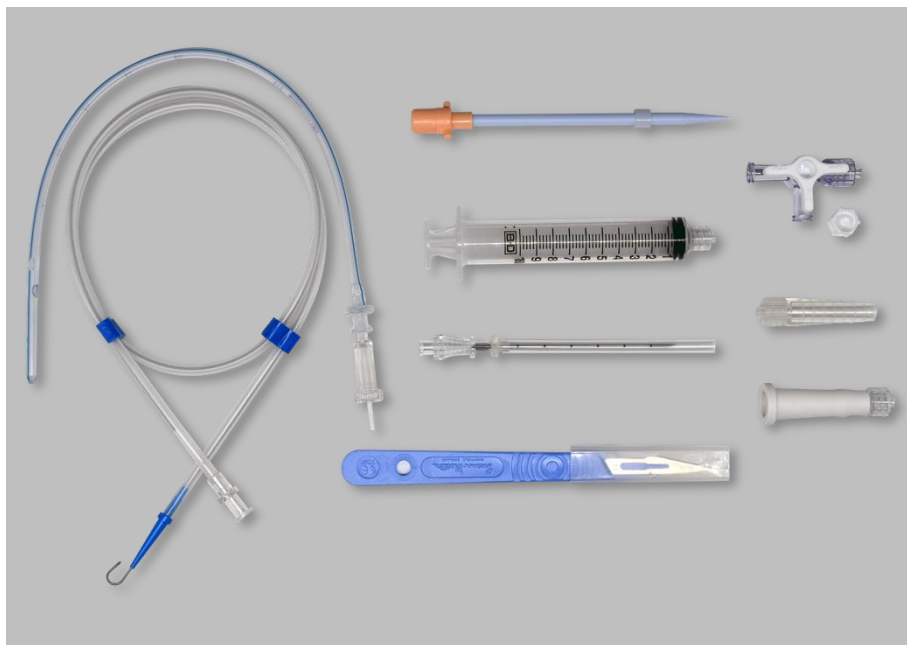

The **Rocket®** Seldinger Chest Drainage Set is indicated for the percutaneous insertion of a pleural drain for fluid and air and is particularly suitable for: ¹

Pneumothorax, spontaneous or iatrogenic
Effusions, benign or malignant, except those with high protein content.

Low viscosity empyema

References:

1. Clinical Evaluation of Seldinger Drainage Kit. Beaumont A C et al Proceedings of Intensive Care Society Meeting. December 1999. 2. Catheters smaller than 24 French gauge can be used for chest drains. Taylor PM. BMJ (1997) 315:186 3. Management of Malignant Pleural Effusions Antunes G, Neville E, Thorax (2000) 55:981-983

The **Rocket®** Seldinger Chest Drainage Set features:

- Improved patient comfort using small bore catheterisation ^{2,3}
- Reliable insertion via Seldinger approach under direct vision, x-ray screening or ultrasound guidance
- Scalpel for initial skin puncture
- Syringe and Seldinger Needle to confirm entry to the pleural cavity and insertion of the guide wire
- Guidewire to maintain the correct track during dilation and drain placement
- Dilator Set to enlarge track through dermal tissue. Each dilator has a movable “Safety Stop” to help prevent over insertion during dilation
- Drainage catheter with radiopaque stripe.
- Three Way Tap & Connector Set including “Stepped” connector and Wide bore luer fittings to permit connection to a wide range of drainage options.

CODE	SIZE	DESCRIPTION
R54544-12-00	12Fg	Scalpel, Seldinger Needle, Guidewire, Dilators, Catheter with radiopaque stripe, syringe, Three way tap, Luer connection set Supplied sterile per carton of 10 units
R54544-18-00	18FG	